

Berkshire



CLEANS, DISINFECTS, AND DEODORIZES.

VersaQuat™

One Step Ready to Use Cleaner & Disinfectant



VersaQuat™ is an EPA-registered, RTU (ready-to-use) one-step disinfectant and cleaner that is an effective solution for cleaning and disinfecting needs. With broad-spectrum, best-in-class efficacy claims, VersaQuat™ delivers consistent quality and intuitive usage, making it an ideal choice for a disinfectant rotation program. In addition, no strong odor, harsh chemicals, or toxins provide flexibility for use and advantages in surface compatibility.


Available in a 32 fl. oz. bottle with trigger sprayer, double bagged to facilitate transfer into the cleanroom.

VersaQuat™ Product Features:

- Kills 99.9% of bacteria, viruses, and fungi
- EPA approved List N - SARS-CoV-2 & variants
- EPA approved List Q - Monkeypox Virus
- Fragrance-free
- Designed for large area cleaning/disinfecting
- Bottles double bagged
- Non-corrosive
- Low odor
- Ready to Use (RTU)
- Made in the USA

VersaQuat™

One Step Ready to Use Cleaner & Disinfectant

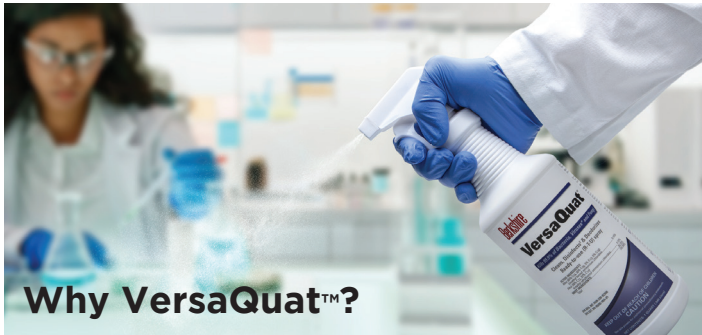

Application Versatility

- Contamination-controlled environments and general institutional/industrial cleaning: Broad-spectrum efficacy (bacteria, viruses, fungi) and one-step cleaner and disinfectant
- Healthcare and long-term care facilities: One step, ready to use effectively and efficiently clean and disinfect nearpatient surfaces and equipment.
- Sensitive surfaces and equipment: Good surface compatibility offers versatility and is not damaging
- Established and trusted chemistry that delivers consistent quality: recommended as part of a disinfectant rotation programs
- Packaging is designed for ease of use in the cleanroom

Quick. Effective. Convenient.

- Non-corrosive
- Non-toxic
- Fragrance-free
- Consistent quality with intuitive usage

Berkshire Product Compatibility

- Can be used with a wide range of Berkshire wipers

Bacteria	Contact Time
<i>Pseudomonas aeruginosa</i> (ATCC 15442)	2 minutes
<i>Staphylococcus aureus</i> (ATCC 6538)	2 minutes
<i>Salmonella enterica</i> (ATCC 10708)	2 minutes
Methicillin Resistant <i>Staphylococcus aureus</i> (MRSA) (ATCC 33592)	2 minutes
Vancomycin Resistant <i>Enterococcus faecium</i> (ATCC 51559)	2 minutes
Multi-Drug Resistant (MDR) <i>Acinetobacter baumannii</i> Ceftazidime and Gentamicin Resistant (ATCC BAA-1605)	2 minutes
Extended Spectrum Beta-Lactamase (ESBL) <i>Escherichia coli</i> (ATCC BAA-196)	2 minutes
Viruses	Contact Time
Avian Influenza Virus (H5N1) NIBRG-14	30 Seconds
Severe Acute Respiratory Syndrome Related Coronavirus 2 USA-WA1/2020 BEI NR-52281	30 Seconds
Human Rotavirus (ATCC VR-2018)	1 minute
Norovirus (Norwalk Virus)	2 minutes
Canine Parvovirus	2 minutes
Rhinovirus Type 14 (ATCC VR-284)	2 minutes
Mildew	Contact Time
<i>Aspergillus niger</i> (ATCC 6275)	10 minutes

* Please refer to product label for comprehensive directions for use.

Quick. Effective. Convenient.

WHERE TO USE:

Hard, non-porous surfaces

DIRECTIONS FOR USE:

GENERAL CLEANING:

Spray soiled area then wipe with a dry paper towel or lint-free cloth.

* DISINFECTION, DEODORIZING:

Apply product with a refillable spray bottle. Hold sprayer six to eight inches from the surface to be treated. Spray area until it is covered with the product. Allow surface to remain visibly wet for 2 minutes. No scrubbing necessary. Allow treated surfaces to air dry or wipe off with a clean cloth, mop or sponge.

* VIRUCIDAL ACTIVITY:

This product is an effective virucide on hard non-porous surfaces against Avian Influenza Virus (H5N1) NIBRG-14, and SARS-Related Coronavirus 2 (SARS-CoV-2) when the treated surface is allowed to remain visibly wet for **30 seconds**. Allow treated surfaces to air dry or wipe off with a clean cloth, mop or sponge.

MILDEWSTAT: To control mold and mildew (*Aspergillus niger*) on pre-cleaned, hard, non-porous surfaces spray surface to be treated making sure to wet completely. Let air dry. Repeat application at weekly intervals or when mildew growth appears.

INGREDIENTS:

ACTIVE INGREDIENT

Alkyl (60% C14, 30% C16, 5% C12, 5% C18)
Dimethyl benzyl ammonium chlorides
Alkyl (68% C12, 32% C14)
dimethyl ethylbenzyl ammonium chlorides 0.15%

EPA REGISTERED

EPA REG. NO: 95747-1-92998

PACKAGING INFORMATION:

UNIT SIZE

32 Fl. Oz.

UNIT OF SALE

6x32 Fl. Oz.

SHELF LIFE 2-YEARS

For order information or to learn more about VersaQuat™, please reach out to one of our Technical Account Managers or visit our website.